

1. Device Name:

The device brand name is Lineage™ Rx PTCA Balloon Dilatation Catheter. The generic name is Rapid Exchange PTCA Balloon Dilatation Catheter.

2. Device Description:

Lineage™ Rx PTCA Balloon Dilatation Catheter which is designed for dilating stenotic atherosclerotic lesions in coronary arteries or bypass grafts and for post-delivery expansion of Balloon expandable stents.

Lineage™ Rx PTCA Balloon Dilatation Catheter consists of a Balloon (dilatation element) near the distal tip with a single radiopaque marker band at the middle of Balloon working length for Balloon sizes 1.25 mm to 1.50 mm and double radiopaque marker bands for balloon sizes 2.00 mm to 4.50 mm, which represent the approximate balloon working length at nominal pressure.

The catheter has a soft tip. The two coaxial lumens permit movement of guide wire and balloon inflation. The two markers on the proximal shaft approximately indicate the exit of the balloon catheter tip from the guiding catheter. (Brachial 90 cm, Femoral 100 cm). The proximal portion of the shaft has PTFE Coating while distal portion of the shaft has hydrophilic coating applied on it, in order to optimize proximal pushability with a smooth transition to a distal shaft trackability.

The compliance chart on the inner and outer label indicates how the balloon diameter increases with increasing pressure. The compliance data is based on in vitro testing of balloons at 37°C temperature.

Table 2.1: Device Component Description

Parameters	Lineage™ Rx PTCA Balloon Dilatation Catheter
Available Balloon Lengths (mm)	6, 9, 12, 14, 15, 17, 20, 25, 30, 33, 38, 41
Available Balloon Diameters (mm)	1.25, 1.50, 2.00, 2.25, 2.50, 2.75, 3.00, 3.25, 3.50, 4.00, 4.50
Delivery system usable length	142 cm
Balloon	Semi Compliant Nylon Balloon. Its Lengths and/or location is defined by single or double swaged radiopaque markers.
Shaft Diameters	Proximal 2.13F Distal 2.7F
Balloon Inflation Pressure	Nominal Pressure: 7 atm. Rated Burst Pressure: 16 atm *RBP 14 atm for balloon dia 4.5 mm.
MGCID	Min Guide Catheter I.D 0.056" / 1.42 mm Guide Catheter Compatibility 5F
Maximum Guide Wire outer Diameter	0.014" / 0.36 mm diameter
Coating	Hydrophilic

3. Indications:

Lineage™ Rx PTCA Balloon Dilatation Catheter is indicated for balloon dilatation of the stenotic portion of a coronary artery or bypass graft stenosis for the purpose of improving myocardial perfusion. Lineage™ Rx PTCA Balloon Dilatation Catheter is also indicated for post-delivery expansion of balloon expandable stents.

4. Contraindications:

- Unprotected left main coronary artery lesions.
- Coronary artery spasm in the absence of a significant stenosis.
- Patients with hypersensitivity or allergies to aspirin, heparin, clopidogrel, ticlopidine or materials such as Pebax, nylon, stainless steel, poly tetra fluoro ethylene, sodium hyaluronate or any contrast media.

5. Warnings:

- For single use only. Do not reuse, reprocess or resterilize. Reuse, reprocessing or re-sterilization may compromise the structural integrity of the device and/or lead to device failure which, in turn, may result in patient injury, illness or death.
- Reuse, Reprocessing or re-sterilization may also create a risk of contamination of the device and/ or cause patient infection or cross-infection, including, but not limited to the transmission of infectious disease (s) from one patient to another. Contamination of the device may lead to injury, illness or death of the patient.
- PTCA should only be performed at hospitals where emergency coronary artery bypass graft surgery can be quickly performed in the event of a potentially injurious or life-threatening complication.
- PTCA in patients who are not acceptable candidates for coronary artery bypass graft surgery requires careful consideration, including potential haemodynamic support during the procedure, as treatment of this patient population carries special risk.

- Use the product before the "Use By" date specified on the product label.
- Guiding catheters used must have lumen sizes that are suitable to accommodate the introduction of Lineage™ Rx PTCA Balloon Dilatation Catheter.

- Do not use if the inner package is open or damaged. Carefully remove the PTCA catheter from the pouch to prevent damage and premature removal of the balloon cover.

- Use only the appropriate Balloon inflation media. Do not use air/any other gaseous media/oil based inflation media/alcohol/organic solvents to inflate the PTCA catheter as this may cause uneven expansion/catheter leakage/loss of lubrication.

- When the catheter is exposed to the vascular system, it should be manipulated while under high-quality fluoroscopic observation. Do not advance or retract catheter unless the balloon is fully deflated under vacuum. If resistance is felt during manipulation, determine the cause of the resistance before proceeding.

- Balloon pressure should not exceed the rated burst pressure. The rated burst pressure is based on the results of in vitro testing. At least 99.9 percent of the balloons, (with a 95 percent confidence) will not burst at or below their rated burst pressure. Use of a pressure monitoring device is recommended to prevent over pressurization.

- To reduce the potential for vessel damage, the inflated diameter of the balloon should approximate the diameter of the vessel just proximal and distal to the stenosis.

- Should any unusual resistance be felt at any time during either lesion access or removal, the entire system must be removed as a single unit.

6. Precautions:

- Prior to angioplasty, examine the PTCA catheter to verify functionality. Ensure that its size and shape are suitable for the specific procedure for which it is to be used.

- Only the physicians trained in the performance of PTCA should use the catheter system.
- Appropriate anticoagulant/antiplatelet therapy should be used during the procedure.
- Prior to reinsertion or withdrawal of the PTCA catheter, wipe the guide wire with saline-soaked gauze to remove blood or other residues.
- After the procedure, anticoagulant therapy should be continued as recommended by the physician.

7. Adverse Effects:

Potential adverse events, which may be associated with the use of a PTCA Balloon Dilatation Catheter, include but are not limited to:

- Acute myocardial infarction.
- Allergic reactions to anti-coagulant and /or antithrombotic therapy/ contrast medium.
- Arrhythmia, including ventricular fibrillation (VF).
- Arteriovenous fistula.
- Coronary embolism.
- Coronary artery spasm.
- Coronary vessel dissection/injury/perforation/rupture.
- Death.
- Emergency or non-emergent Coronary Artery Bypass Graft Surgery.
- Hematoma.
- Hemorrhage, requiring transfusion.
- Hypotension / Hypertension.
- Infection and / or pain at the access site.
- Restenosis of treated segment.
- Total occlusion of coronary artery/bypass graft.
- Unstable angina pectoris.
- Stroke, air embolism and embolization or fragmentation of thrombotic or atherosclerotic material

8. How Supplied:

Sterile: This device is sterilized with Ethylene oxide, Non pyrogenic.
Do not use if the package is open or damaged.

Contents: One (1) Lineage™ - Rx PTCA Balloon Dilatation Catheter housed in a protective circular hoop dispenser,

One (1) Instructions for use, One (1) Re - grooming sheath for refolding balloon, if needed, One (1) Flushing needle, One (1) Looper clip to hold the catheter in a coil configuration, One (1) Compliance chart.


See sections 9.5 and 9.6 for the Re-grooming sheath and the Looper Clip instructions for use.

See the figure below for the placement of the contents in Lineage™ package.

Storage: Store in a cool, dry and dark place in its original packaging.

Note:

- A - Lineage™ - in Hoop Dispenser
- B - Flushing Needle
- C - Re-grooming Sheath
- D - Looper Clip
- E - Compliance chart

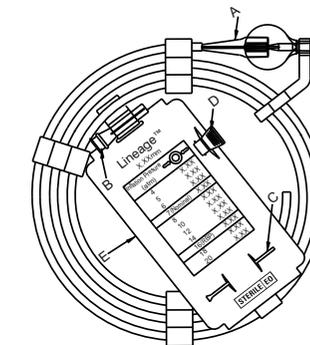

9. Operator's Manual:

The following instructions provide guidance but do not obviate the necessity of training in the use of the device.

- The PTCA procedure should be carried out according to standard PTCA guidelines.
- Read Instructions for Use.

9.1 Preparation:

- Open product box to remove the inner Tyvek pouch containing a hoop dispenser with sterile Lineage™ Rx PTCA Balloon Dilatation Catheter.

Note: Do not use if the package is open or damaged. Remove carefully to avoid damage.

- Remove protective sheath and stylet from the tip.
- Prior to using Lineage™ Rx PTCA Balloon Dilatation Catheter, carefully inspect for bends, kinks, and other damages. Ensure that the radiopaque marker(s) are placed in the balloon.
- Insert the flushing needle into the distal tip of the PTCA Catheter. Flush the guide wire lumen of the Catheter with sterile heparinised solution by gently applying pressure with the syringe for at least 10 seconds. The saline should exit from the guide wire port. Remove the syringe while flushing the needle.
- Attach a stopcock to the Catheter's inflation hub.
- Attach an appropriate size of syringe filled with a minimum of 3cc of 50% solution of contrast medium in sterile heparinised saline to the stopcock.
- Open the stopcock and induce negative pressure, pulling the syringe plunger as far back as possible without dislodging it from the syringe barrel. Aspirate till the air is fully evacuated.
- Gently release the negative pressure to fill the Catheter lumen with the solution.
- Close the stopcock, remove the syringe and purge the air.
- Repeat steps (g) through (i) if, necessary till there are no air bubbles left.
- Fill an Angioplasty Inflation Device with a minimum of 3cc of a 50% solution of contrast medium in sterile heparinised saline and attach it to the purged Catheter.

Note: Purge the Angioplasty Inflation Device with air prior to connecting to the Catheter.

- Induce negative pressure and close the stopcock. Purge any air from the Angioplasty Inflation Device through the stopcock. Leave the Catheter under negative pressure until it is ready to be used.

9.2. Assembly and Insertion Procedure for PTCA:

- Attach a haemostatic valve connector to the proximal luer of the guide catheter already positioned within the vasculature.
- Insert the guide wire carefully through the haemostatic valve connector. While observing under fluoroscopy, cross the lesion with the guide wire according to the accepted PTCA techniques.
- Gradually tighten the haemostatic valve to avoid back flow of blood.
- Backload the distal end of the guide wire into the tip of the Balloon Catheter. The guide wire will exit through the guide wire exit port.

Note: When back loading the Catheter, provide adequate support to all shaft segments. Do not try to straighten a kinked Catheter.

- Open the valve to facilitate easy movement of the Balloon Catheter over the guide wire.
- Advance the PTCA Catheter until appropriate shaft markers correspond to the haemostatic valve hub. This indicates that the tip of the balloon catheter has reached the distal end of the guiding Catheter.

Warning: If a Tuohy-Borst type adjustable haemostatic valve is used, avoid over tightening since this may restrict the flow of contrast medium in and out of the inflation lumen, thereby slowing inflation/deflation.

Caution: To avoid shaft movement restrictions, adjust haemostatic valve when advancing/withdrawing Catheter, allowing for changes in the shaft diameter.

Caution: During the advance, ensure that negative pressure is maintained in the Inflation Device to keep the balloon fully deflated.

Caution: Do not advance or withdraw the PTCA Catheter within the coronary vasculature unless it is preceded by a guide wire.

Caution: Care should be taken when inserting the PTCA Catheter into the haemostatic valve in order to avoid kinking.

- Determine the correct position of the Balloon Catheter under fluoroscopy using radiopaque markers.
- If resistance is encountered during advance, do not force the passage. Determine the cause of resistance and take corrective actions.

9.3 Inflation:

- Dilate the stenotic lesion following standard PTCA techniques.

Caution: Use of oversized balloon may cause dissection. Do not exceed the Rated Burst Pressure mentioned on the product labels.

- If needed, carry out repeat inflations. Maintain negative pressure between inflations.
- After each subsequent inflation, determine the distal blood flow.
- If significant stenosis persists, the lesion may require successive inflations/further treatment.
- Determine the results under fluoroscopy.

9.4 Catheter Exchange/Withdrawal and Disassembly Procedure:

- Apply negative pressure to the inflation device to achieve complete deflation of the Balloon. Confirm this under fluoroscopy.
- Loosen the haemostatic valve to facilitate removal.
- While holding the guide wire and the haemostatic valve with one hand, gradually remove the balloon catheter from the guide catheter.

- In case of exchange procedure, follow steps (e) through (i). For withdrawal or disassembly, follow steps (j) and (k).
- Maintain the guide wire position if reinsertion/further procedure is planned.
- Withdraw the PTCA Catheter till the guide wire exit port is exposed.
- Slide the remainder of the PTCA catheter along the guide wire until the tip of the PTCA catheter exits the haemostatic valve. Retighten the haemostatic valve onto the guide wire to hold it securely in place.
- Completely remove the PTCA Catheter from the guide wire. Clean with gauze soaked in sterile normal saline to wipe residues of blood and other substances. The mandrill and the protective sheath can be used to cover the balloon by sliding the distal pin into the tip and the protective sheath over the balloon with the flared side of the protector sheath first.
- Prepare and insert the next catheter according to manufacturer's instructions. If reinserting the same Lineage™ Rx PTCA Balloon Dilatation Catheter, flush the guide wire lumen of the Lineage™ Rx PTCA Balloon Dilatation Catheter using the flushing needle (as described in step (d) under 9.1 Preparation). Before use, examine the PTCA Catheter carefully for any damage. If damaged, do not reuse.
- Withdraw the deflated PTCA Catheter and guide wire from the guiding catheter.
- Remove the guiding catheter from the vasculature following an accepted technique for removal.

9.5 Instruction for use of Accessories:

9.5.1 Re-grooming Sheath:

The Re-grooming sheath is an accessory component that allows the balloon to be re-folded if needed during the procedure.

Instruction for use:

- Deflate the balloon by applying negative pressure to the inflation device and maintain under vacuum.
- Visually inspect the balloon to confirm that it is fully deflated.
- Carefully load the stylet back through the distal tip of the catheter then load the flared end of the re-grooming sheath over the stylet.
- Hold the catheter just proximal to the balloon and push the re-grooming sheath over the balloon in a gentle twisting motion until the entire balloon is covered.
- Gently remove the re-grooming sheath and the stylet from the balloon catheter.
- Visually inspect the balloon for any potential damage. Discard the balloon catheter if there is any visible damage on the balloon.

9.5.2 Flushing needle


This needle is used to flush the guide wire lumen.

9.5.3 Looper clip

The Looper clip is an accessory component that allows the catheter to be fastened in a coiled configuration for ease of handling during use (See the figure below).

Instruction for use:

- Form the catheter into a single or double loop as needed.
- Press down the arms of the Looper clip to open the clipping section for retention of the hypotube.
- Release the arms of the Looper clip with the hypotube inside to secure the catheter in a coiled configuration. Discard the catheter if any bend or kink is observed upon placement or removal of the Looper clip.

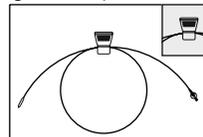

9.5.4 Compliance Chart:

To provide balloon compliance information.

9.6 Disposal:

After usage, dispose of the balloon catheter and accessories in accordance with accepted regulations for medical waste management.

Note: The literature should be consulted for updates to the current PTCA procedural guidelines.

10. Balloon Compliance Data:

Table 10.1: Balloon compliance chart

Pressure (atm)	1.25 mm	1.50 mm	2.00 mm	2.25 mm	2.50 mm	2.75 mm	3.00 mm	3.25 mm	3.50 mm	4.00 mm	4.50 mm
4	1.16	1.40	1.90	2.11	2.37	2.59	2.85	3.09	3.36	3.88	4.30
5	1.19	1.44	1.95	2.15	2.43	2.64	2.91	3.15	3.41	3.94	4.35
6	1.21	1.48	1.98	2.19	2.46	2.69	2.96	3.20	3.46	3.98	4.42
7	1.25	1.50	2.00	2.25	2.50	2.75	3.00	3.25	3.50	4.00	4.50
8	1.26	1.53	2.04	2.29	2.53	2.77	3.03	3.30	3.53	4.05	4.54
10	1.30	1.58	2.08	2.35	2.58	2.82	3.08	3.37	3.59	4.11	4.60
12	1.33	1.63	2.12	2.39	2.63	2.87	3.12	3.43	3.63	4.15	4.66
14*	1.36	1.67	2.16	2.41	2.67	2.91	3.16	3.49	3.67	4.20	4.72
16	1.38	1.72	2.20	2.43	2.71	2.95	3.20	3.54	3.72	4.25	

Nominal Dilatation Pressure 7 atm for all diameters. **Grey background:** Nominal pressure, **Black background:** RBP (Rated Burst Pressure) *RBP 14 atm for Balloon diameter 4.5mm.


11. Disclaimer of Warranty and Limitation of Remedy:

There is no express or implied warranty, including without limitation any implied warranty of merchantability of fitness for a particular purpose, on the Meril Life Sciences Pvt. Ltd. product(s) described in this publication. Under no circumstances shall Meril Life Sciences Pvt. Ltd. be liable for any direct, indirect, incidental or consequential damages. Some jurisdictions do not allow the exclusion of or limitations on an implied warranty. For jurisdictions which do not allow the exclusion or limitations of incidental or consequential damages, some of the above exclusions may not apply. The patient may also have other rights, which vary from jurisdiction to jurisdiction.

Descriptions or specifications in Meril Life Sciences Pvt. Ltd. printed matter, including this publications, are meant solely to generally describe the product at the time of manufacture and do not constitute any express warranties.

Meril Life Sciences Pvt. Ltd. will not be responsible for any direct, indirect, incidental, or consequential damages resulting from reuse of the product.

Symbols used in labeling

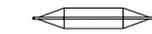

Balloon outer diameter



Contains one unit



Do not use if package is damaged



Reference number



Max. guide wire diameter



Consult instructions for use




Sterilized using ethylene oxide



Do not re-use



Authorized representative in the European community

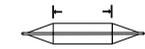

Balloon length



Keep dry



Non-Pyrogenic



Lot number



Min guide catheter I.D.



Use-by date



Manufacturer



Date of manufacture



Do not re-sterilize

Meril

Manufactured by:
Meril Life Sciences Pvt. Ltd.
Muktanand Marg,
Chala, Vapi 396191,
Gujarat, India.
Web: www.merillife.com



Obelis s.a.
Bd., General Wahis 53,
1030, Brussels, Belgium.
Tel.: +32.2.732.5954
Fax: +32.2.732.6003
E-mail: mail@obelis.net



Customer Care Contact:
Tel.: +91 (260) 240 8000
E-mail: askinfo@merillife.com

An ISO 13485 & GMP Certified Company